




Battery Jigs & Holders

Battery & Coin Cell Jig

- Holds different sizes of cylindrical cell, coin cell, pouch cell, prismatic cell
- Wide contact point with noble coated contact area
- 4 contact point type (Kelvin probe) is available to minimize voltage drop for high current application.
- Individual channel operation is available.
- Can be moved to fit cell size

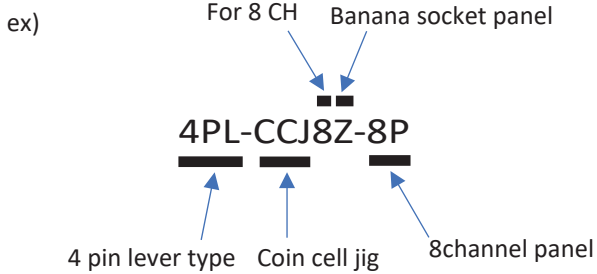

* Ordering Guide

	Part No.	Description
1st	4PL	4 Pin Lever type
	4PK	4 Pin Knob type
2nd	CCJ	Coin cell jig
	UCJM	Universal cell jig
	UCJH	Universal jig for high current
	PCJ	Pouch cell jig
3rd	Channel No.	Channel Quantity
4th		Normal type
	Z	Banana socket panel
5th	4P	4ch per panel
	8P	8ch per panel

4PL-CCJ8Z-8P



- 8 channel Coin Cell Jig
- 4pin lever type
- banana socket panel (without jig cable)
- 8ch per pannel

4PL-UCJM8Z-8P



- 8 channel Universal Cell Jig
- 4 pin lever type
- banana socket panel (without jig cable)
- 8ch per pannel

4PK-UCJH4-4P



- 4 channel universal jig for high current
- 4 pin knob type
- 4ch per panel

4PL-PCJ8-4P



- 8 channel Pouch Cell Jig
- 4 pin lever type
- 4ch per pannel
- including 3m jig cables

UCJ1A



- 1 channel universal jig
- 4pin lever type

PCJ1



- 1 channel Pouch Cell Jig
- 4pin lever type
- For small size pouch cell

PRCJ1



- 1 channel Prismatic Cell Jig
- 4 pin probe knob
- For prismatic battery or pouch cell

High current Cylindrical Battery Holder

- *These Battery holders are for single cell high current cylindrical battery test*
- *Kelvin type 4 probe connection to minimize voltage drop by cable resistance and/or contact resistance etc.*
- *Maximum current : 30A*
- *Maximum battery diameter : 30 mm*
- *Minimum battery contact diameter : 14mm*
- *Maximum battery length : 65 mm or 100 mm*
- *Current path diameter : 14mm*
- *Minimum battery length : 1 mm*
- *Cable Connectors : 4ea of 4mm banana*
- *Size*
 - HCCBJ65L : 136x24x43 mm, 208g*
 - HCCBJ100 : 172x24x43 mm, 240g*
- *Part No.*
 - max. battery length 65mm : HCCBJ65L*
 - max. battery length 100mm : HCCBJ100*



HCCBJ65L



HCCBJ100L

Coin Cell Holder

For WMPG/WBCS Series

- *Direct connect to cell connector*
- *Part No.*
 - for standard current model : CCH2*
 - WMPG1000S, WBCS3000S*
 - for low current model : CCH2L*
 - WMPG1000Ls/Le, WBCS3000Ls/Le/Lx*



For ZIVE Series

- *D-SUB connector type*
- *Part No.*
 - for 20mm dia. coin cell : CCH3-20*
 - ZIVE SP1, SP2, SP3, MP1, MP2, BP2*
 - for 24mm dia. coin cell : CCH3-24*
 - ZIVE SP1, SP2, SP3, MP1, MP2, BP2*



• • • **Related Instruments**

Battery Test System



Impedance Analyzer



Battery Impedance Analyzer



WonATech Co., Ltd.
 7, Neunganmal 1-gil, Seocho-gu,
 Seoul, 06801, Korea
 Phone: +82-2-578-6516
 Fax: +82-576-2635
 E-mail: sales@wonatech.com
 Website: www.wonatech.com

Local Distributor

